


INVACARE®  
**PLATINUM™ 10**  
OXYGEN CONCENTRATOR

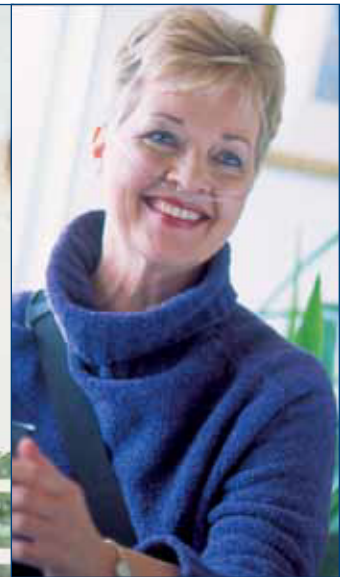

---

The Invacare® Platinum™ 10 Oxygen Concentrator provides an attractive alternative to liquid systems and is known for reliability, performance and value. The 10 liter unit features quiet operation, a top handle and bottom recesses for easy handling and a side filter access door. In addition, the concentrator is compatible with the Invacare® HomeFill® Oxygen System and is available with an Invacare® SensO<sub>2</sub>® Oxygen Purity Sensor.

---

## Features

Pressure-based system is superior to timing-based systems  
 Self diagnostic electronics minimize troubleshooting time  
 Top handle for ease of mobility  
 Battery-free power loss alarm, Invacare® Sieve-GARD™ Monitor protects sieve from degradation, compressor 35 psi pressure-relief valve and low flow alarm

## Provider Benefits

Provides continuous flow up to 10 liters per minute as an alternative to liquid systems  
 All Platinum 10 stationary concentrators include a HomeFill system port for inventory flexibility  
 Patient-accessible reset button helps minimize emergency service calls

## Patient Benefits

User friendly design with easy to see control panel, humidifier bottle holder and filter access  
 HomeFill system compatibility allows patients with flow rates up to 5LPM to enjoy life's experiences with small, refillable cylinders  
 Reliable concentrator provides unlimited, worry-free oxygen supply

## HCPCS Code: E1390



1. Invacare® HomeFill® Oxygen System connection
2. Easy to maintain filter
3. Also available with Invacare® SensO2® Oxygen Purity Sensor



**Invacare Corporation**  
 www.invacare.com

### USA

One Invacare Way  
 Elyria, Ohio  
 44035-4190  
 (800) 333-6900

### Canada

570 Matheson Blvd. E.,  
 Unit 8  
 Mississauga, Ontario  
 L4Z 4G4 Canada  
 (800) 668-5324

©2013

Invacare Corporation.  
 All rights reserved.  
 Trademarks are identified by the symbols ™ and ®. All trademarks are owned by or licensed to Invacare Corporation unless otherwise noted. Specifications are subject to change without notification.  
 Form No. 11-190  
 Rev. 04/13  
 130376



**Yes, you can.®**

Model #	Dimensions	Weight	Shipping Weight	Quantity	Oxygen Flow	Flow Range	Maximum Outlet Pressure	Pressure Relief Mechanism
IRC10LX	26 3/8" H	54 lb.	60 lb.	1 Per Case	10 LPM	2-10 LPM	9+/- 0.5 PSI	35 PSI
IRC10LXO2	18 3/8" W 14 3/8" D			18 Per Pallet				+/- 3.5 PSI

O <sub>2</sub>	Filter	Electrical Requirements	Rated Current Input	Average Power Consumption	Sound Level	Limited Warranty
Oxygen Purity	Cabinet	120 VAC	4.0 A	585 W	50 db. Average	3 Years
94%-87%	At All Flow Rates	Compressor Inlet, +/- 10%, 60 Hz Outlet HEPA (132VAC/108VAC)				

All references to HCPCS codes: Providers are responsible for determining the appropriate billing codes when submitting claims to the Medicare program, and should consult an attorney or other advisor to discuss specific situations in further detail.